



BEZNOSKA

We bring back joy to movement



Cementless Revision Stem Type SF



Revision Implants of Hip Joint

ARTHROPLASTY


■ Introduction

The design of the cementless revision stem is based on the type S.F. total hip replacement stem, which is exceptionally also used for re-implantation, even though this implant is intended for primo-implantation. This concept allows the use of the majority of instruments from the instrumentarium intended for type S.F. stem insertion. The designers made use of their experience with type S.F. cementless stems, and the manufacturing process has utilized the most sophisticated technologies.

When using instruments supplied by the manufacturer, straightforward implantation and perfect fixation of the implants can be assured. The stem assortment enables the majority of cases that may arise in revision surgery of total hip replacement to be solved.

This publication should serve as surgical guidance for the particular implant and instrumentarium. For the sake of conciseness, it is focused only on the implantation of the given type of endoprosthesis, and it is assumed that the surgeon and other theatre staff are perfectly acquainted with general surgical principles of hip joint arthroplasty.

The aim of this publication is to provide surgeons and scrub nurses with a quick reference guide on the correct use of individual components of the instrumentarium, so that optimum outcomes are achieved and – last, but not the least – unnecessary damage or degradation of the instrumentarium and/or the implant are avoided. By no means is this intended as a manual for surgical technique.

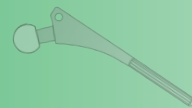


■ Femoral Stem

The hip joint revision stem intended for implantation without bone cement is manufactured from Ti6Al4V titanium alloy. It is supplied in two lengths – 190 mm and 230 mm, and in six diameters of the cylindrical groove stem – 13, 15, 17, 19, 21 and 23 mm.

The stem has 9 – 12 grooves (depending on the stem diameter), which ensure the rotatory stability of the stem.

The upper part of the stem is coated with plasma-sprayed titanium in the length of 95 - 115 mm. The lower part of the stem is cut in three directions to reduce its stiffness. The endoprosthesis neck is ended with a 12/14 cone (EURO), and the CD angle of the endoprosthesis is 135 grades.

Type SF Cementless Revision Stem – 12/14 Conical Neck


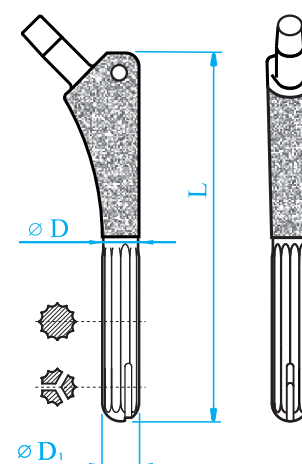



Straight Stem, L = 190		
Size	Ø D / Ø D ₁ (mm)	Order Number
0	10/11	322050
1	12/13	322051
2	14/15	322052
3	16/17	322053
4	18/19	322054
5	20/21	322055
6	22/23	322056

Straight Stem, L = 230		
Size	Ø D / Ø D ₁ (mm)	Order Number
0	10/11	322060
1	12/13	322061
2	14/15	322062
3	16/17	322063
4	18/19	322064
5	20/21	322065
6	22/23	322066

■ Instruments

The instrumentarium is placed in two trays, enabling a well-arranged layout for the instruments during surgery, and for their preparation, sterilisation and storage. The layout of the instruments in trays is pictured in the following figures. During transportation, the trays are placed in a Styrofoam container and then placed in a cardboard outer package.

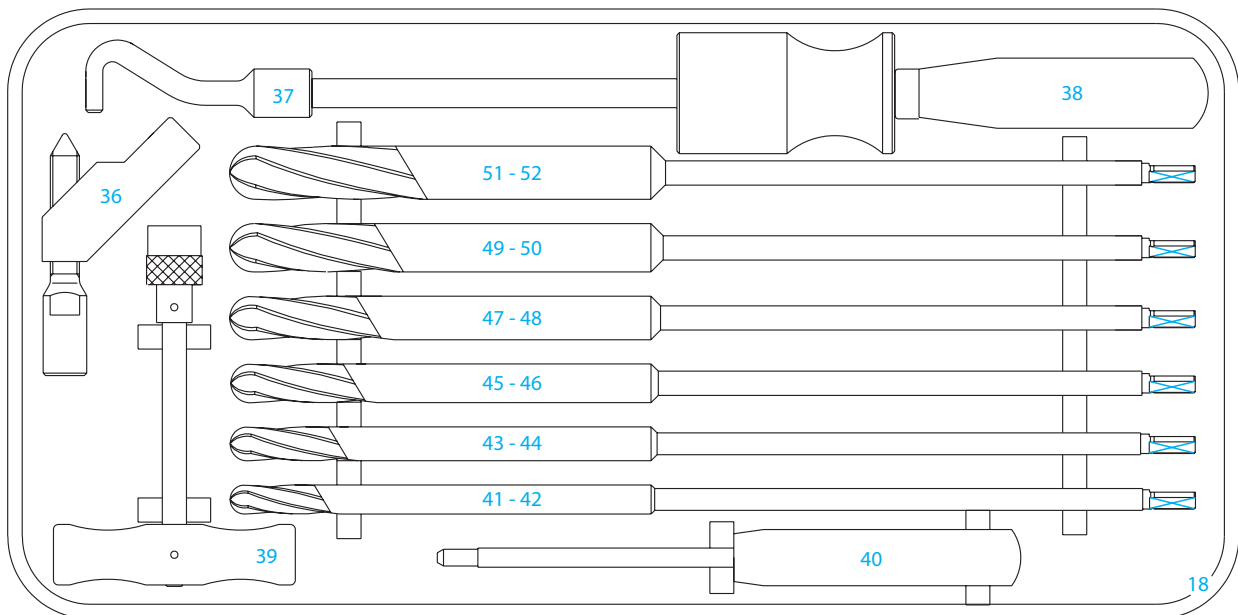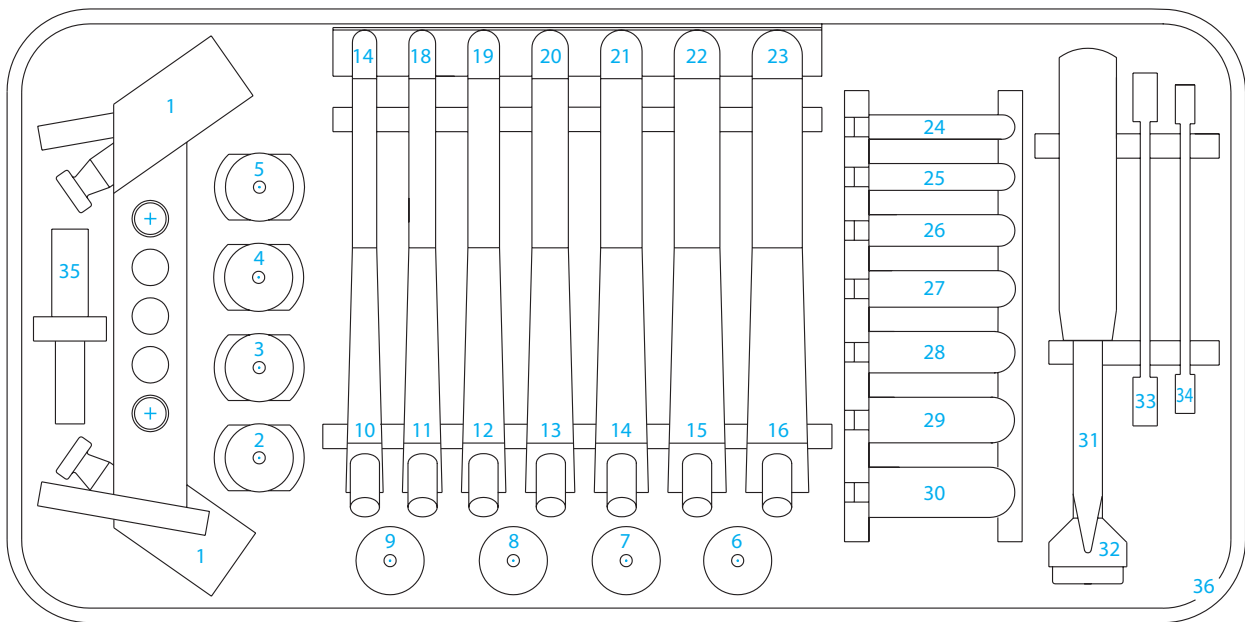

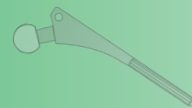

Surgical Technique

General Principles for Implantation of Type S.F. Cementless Revision Stem

Prior to any revision surgery of total hip arthroplasty, a decision must be made whether this type of implant is suitable for the given case. Then it is necessary to carry out preoperative planning, allowing the size of the femoral component to be determined. Preoperative planning must include a pelvic x-ray, with both hips in A/P views, and axial views of both hip joints to determine the width of the metaphysis. The planning of the stem size is performed using transparent templates provided by the manufacturer. Place these templates on the x-ray images of the same scale and determine the length and diameter of the revision endoprosthesis stem.

Perform the revision surgery using an anterolateral approach by Bauer. During dissection around the upper end of the femur, it is convenient to use an electric scalpel.

After removal of the original implant and subsequently any residual bone cement, it is necessary to send any present effusion, soft tissue samples, and the removed implant for microbiological analysis and testing for antibiotic sensitivity; the implant has to be tested to determine microbiological findings in the film covering the implant surface.

After the bone marrow canal is prepared according to the mentioned surgical technique, insert the revision stem - in case of no or minimum damage to the proximal femur - directly to the bed prepared by the instruments. In case of more extensive damage, spongy debris or bone grafts should be used.


Primary fixation of the stem is ensured by insertion with a 0.5 mm overlap. In the proximal part of the stem, it is an overlap of the porous sprayed titanium, and in a cylindrical stem, it is an overlap of the groove struts, which ensures rotatory stability.

This primary fixation achieved during surgery is enhanced by ingrowing of the bone tissue, which is called secondary fixation. This fixation develops gradually during the period of 10 – 12 weeks after implantation. After this period, the implant is usually completely joined with the bone.

Instruments Needed for Insertion of SF Cementless Revision Stem

Straight Stem, L = 190		
Stem Diameter (mm)	Cylindrical Reamer Diameter (mm)	Rasp Size
10/11	10	0
12/13	12	1
14/15	14	2
16/17	16	3
18/19	18	4
20/21	20	5
22/23	22	6

Straight Stem, L = 230		
Size	Ø D / Ø D ₁ (mm)	Order Number
0	10/11	322060
1	12/13	322061
2	14/15	322062
3	16/17	322063
4	18/19	322064
5	20/21	322065
6	22/23	322066

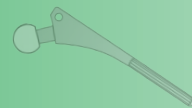

3. Preliminary Articulation Test

Perform preliminary articulation by using test heads ② - ⑤ placed on the cylindrical end of the rasp. The test heads are available with a short, medium, long or extra-long neck (Fig. 3). Once the rasp is removed from the medullary canal, the handle with the ejector pin must be re-mounted to the rasp ③⑤.

4. Insertion of Endoprosthesis

To insert the cementless revision stem in the femur, use a stem impactor ③①. To enable it to be fitted on the endoprosthesis, there is a 4 mm diameter hole in the stem axis. Insert the component so that it sits firmly on the femoral calcar (Fig. 4).

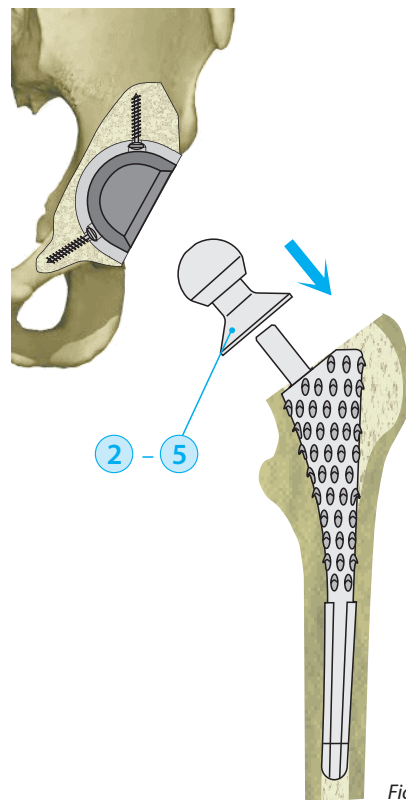


Fig. 3a

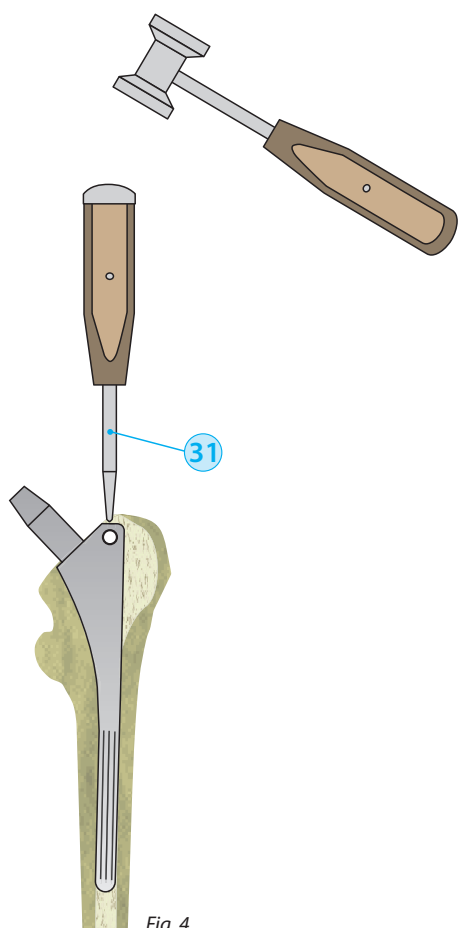


Fig. 4

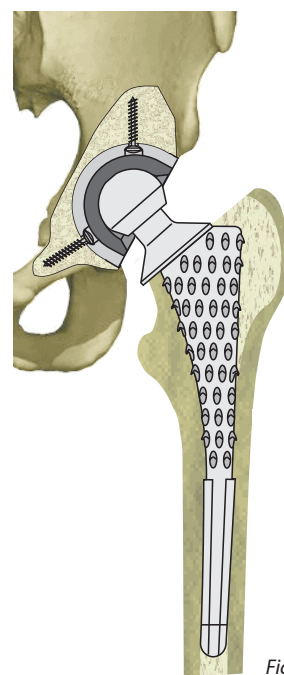


Fig. 3b

5. Final Articulation Test

To carry out the final articulation test, use one of the four plastic test heads – S, M, L, XL ⑥ - ⑨. Based on this test, select the corresponding metallic or ceramic head (Fig. 5).

6. Head Articulation

For articulation of the endoprosthesis head, use a head positioner ③② (Fig. 6).

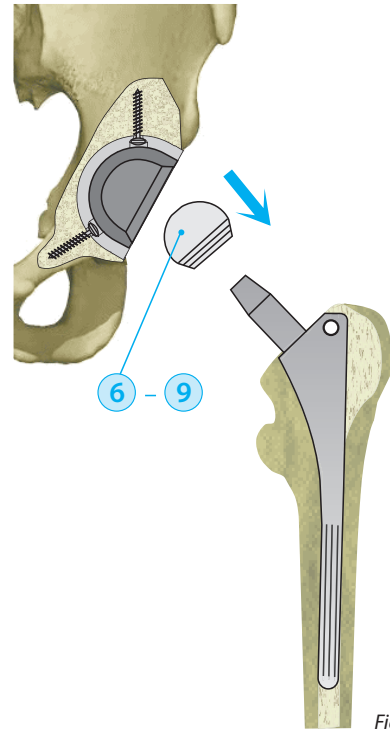


Fig. 5


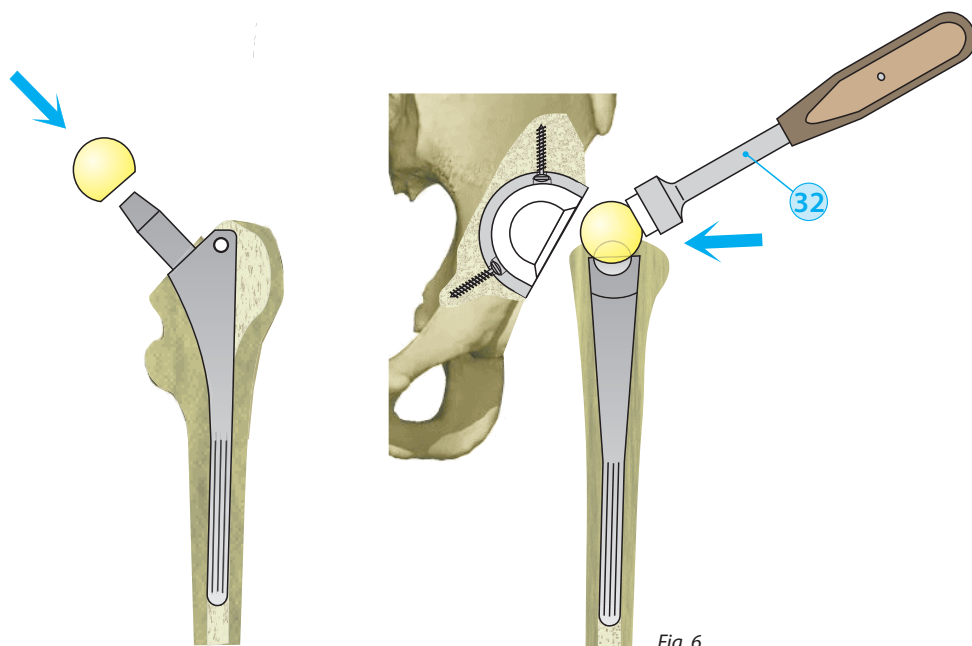


Fig. 6

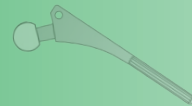

7. Final Result

The final result is visualized in Fig. 7.

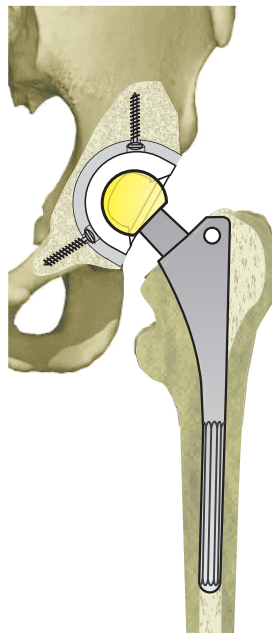


Fig. 7




BEZNOSKA


We bring back joy to movement



■ Contact

Export:

tel.:  +420 312 811 221

 +420 312 811 225

export@beznoska.cz

CE 1014

BEZNOSKA, s.r.o.
Dělnická 2727, 272 01 Kladno, Czech Republic
www.beznoska.com, export@beznoska.cz

tel.: +420 312 660 670
GSM: +420 602 666 503
fax: +420 312 660 216
+420 312 662 464